

MIKROMEISTER®
MASTERING LIQUIDS

MIKROMEISTER HF MASTER OF HYDROFLURIC ACID

- highest Quality made in Germany
- robust for extreme applications
- self priming system
- very high resistance
- wiping seal eliminates swelling and freezing of the piston
- autoclavable 121°C




MIKROMEISTER HF

ROBUST PRECISION

The MIKROMEISTER® HF is the ideal bottle top dispenser for the precise dosing of special aggressive acids and alkalis. Thanks to the use of high-purity materials, it is also suitable for dosing hydrofluoric acid. All wetted parts are made of PTFE, PFA, ETFE and FEP. High-purity ceramic and platinum-iridium are built into the valves.

Scope of supply

MIKROMEISTER HF, 1 filling tube,
3 adapters of PP (GL 28/40/45),
operation manual, certificate of calibration.

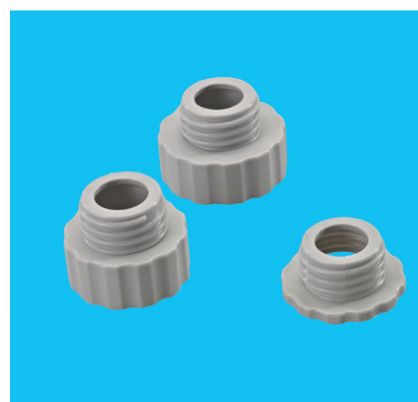

Cat. No. 4207

Volume	Division	Tolerance	Deviation	Thread
ml	ml	+/- %	%	GL
1-10ml	0,2	0,5%	0,1%	28; 40; 45

MIKROMEISTER®
MASTERING LIQUIDS

MIKROMEISTER GmbH
Reichenäcker 11
97877 Wertheim | Germany
Phone +49 9342 8596181
Fax +49 9342 8596183
info@mikromeister.de
www.mikromeister.de